


The Thermo Scientific Orion Star A215 pH/Conductivity Benchtop Meter works for any lab or application. Its durability and dependability make it ideal for every sample, every day.

Thermo Scientific Orion Star A215 pH/Conductivity Benchtop Meter

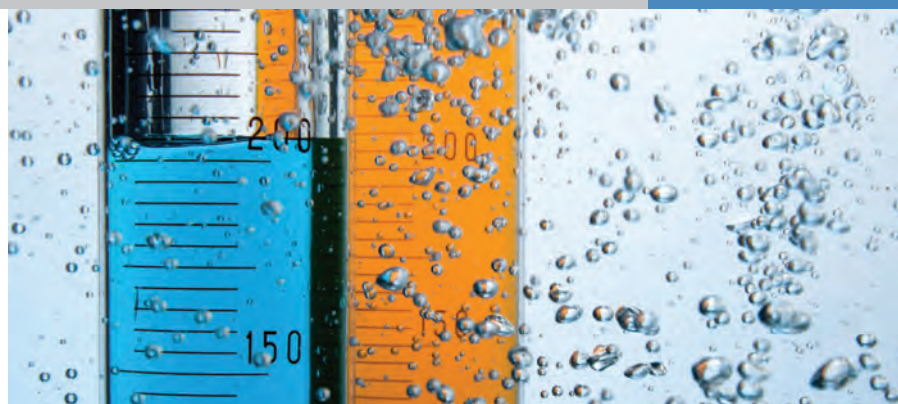

The Thermo Scientific Orion Star A215 pH/Conductivity Benchtop Meter is the ideal choice for any lab. Get the information you need quickly and easily from the large, backlit graphic LCD display. Informative display shows individual pH/mV/RmV or conductivity/TDS/salinity/resistivity results with temperature, or all results at once. Additional helpful information includes electrode status, time, date, sample ID, user ID and calibration. For advanced users, features such as stability and averaging choices provide additional options.

Features and Benefits

- Anyone can use this meter with its plain language screen prompts, soft keys that update for easy selection and a multilanguage interface that supports English, Spanish, French, Italian, German and Chinese
- Don't miss a reading – AUTO-READ™ locks in the stable reading on your screen, ready indicator alerts when readings are stable and timed reading gathers data in specific time intervals
- Up to 5 point pH calibration with automatic recognition for USA/NIST and DIN buffers and calibration editing to fix errors without a complete recalibration
- Selectable cell constant allows for use with 2- or 4-cell conductivity cells
- For accurate conductivity results, selectable reading reference temperatures of 5, 10, 15, 20 or 25 °C with linear, non-linear, nLFu or EP curve options
- Choice of practical salinity or natural sea water curves for salinity readings
- Linear or ISO/EN 27888 curves for TDS
- Non-volatile memory holds up to 2000 data points with time and date stamp
- Easily transfer data and keep meter software up-to-date with the USB and RS232 ports and complimentary data analysis software
- Mix samples without a stir plate with direct control of the Orion Star stirrer probe (096019, sold separately)
- Included electrode arm and newly-designed probe holder make it easier to maintain and place probes into samples
- Works with almost every AC power source with the included universal power adapter, or use four AA batteries (sold separately) to run the meter on DC power
- IP54-rated housing handles splashes and is wall-mountable if table space is a problem
- 3 year meter replacement warranty

Thermo Scientific Orion Star A215 pH/Conductivity Benchtop Meter

Specifications



pH	Range	-2.000 to 20.000
	Resolution	0.1, 0.01, 0.001
	Relative Accuracy	±0.002
	Calibration Points	Up to 5
	Calibration Editing	Yes
mV/RmV	Range	±2000.0 mV
	Resolution	0.1
	Relative Accuracy	±0.2 mV or ±0.05 % of reading whichever is greater
	E_H ORP Mode	Yes
Conductivity	Range	0.001 µS to 3000 mS
	Resolution	0.001 µS minimum; 4 significant figures minimum
	Relative Accuracy	0.5 % reading ±1 digit
	Reference Temperature	5, 10, 15, 20, 25 °C (default)
	Temperature Compensation	Linear (0 to 10.0 %/°C), nLn, nLFu, EP
	Compatible	0.001 to 199.9
	Cell Constants	
	Calibration Points	Up to 5 points
	Calibration Editing	Yes
	Resistivity	Range
Resolution		2 ohms-cm
Relative Accuracy		0.5 % reading ±1 digit
Salinity	Type	Practical salinity or natural sea water
	Range	0.01 to 80.0 ppt NaCl equivalent; 0.01 to 42 ppt practical salinity
	Resolution	0.01
	Relative Accuracy	±0.1
TDS	Range	1 to 200 ppt
	Resolution	4 significant digits
	Relative Accuracy	0.5 % reading ±1 digit
	TDS Factor Range	Linear 0.01 to 10.00, default 0.49
Temperature	Range	-5 to 105 °C, 23 to 221 °F
	Resolution	0.1
	Relative Accuracy	±0.1
	Offset Calibration	1 point
Datalogging	Memory	2000 with time and date stamp
	Log Function	Manual, ready (includes AUTO-READ™), timed
	Log Edit	Delete individual points, group of points or all
Inputs	pH Electrode	BNC, reference pin
	Conductivity with Temperature Probe or ATC Probe	8-pin mini-DIN
Output		RS232, USB
Power	AC Adapter	Included – universal, 100-240 VAC
	Battery Power	Optional – 4 AAs
	Battery Life	800 hrs

Ordering Information

Cat. No.	Description
STARA2150	Orion STAR A215 pH/Conductivity Benchtop Multiparameter Meter - Electrode arm with redesigned holder - Universal power adapter
STARA2155	Orion STAR A215 pH/Conductivity Benchtop Multiparameter Meter Kit - 8157BNUMD Orion ROSS Ultra Triode pH/ATC electrode - 013005MD Orion DuraProbe conductivity cell - ROSS solution kit (475 mL each of pH 4, 7 and 10 buffers; storage solution; cleaning solution; and pH electrode storage bottle) - 011007 Orion 1413 µS conductivity standard, 5x60 mL - Electrode arm with redesigned holder - Universal power adapter
Accessories	
096019	Orion Star stirrer probe
STARA-BEA	Benchtop electrode arm for Orion Star A-series meters
STARA-HB	Freestanding base for use with Orion Star A-series benchtop electrode arm
8302BNUMD	Orion ROSS Triode 3-in-1 pH/ATC probe, refillable, glass body
8102BNUWP	Orion ROSS Ultra refillable pH electrode, glass body
927007MD	Orion ATC probe, stainless steel
810199	Orion ROSS All-in-One pH buffer kit
910104	Orion pH 4.01 buffer, 475 mL
910107	Orion pH 7.00 buffer, 475 mL
910110	Orion pH 4.01 buffer, 475 mL
013005MD	Orion DuraProbe conductivity cell, 4-cell, K=0.475
013016MD	Orion 2-electrode conductivity cell with flow cell, K=0.1
011008	Orion conductivity standard, 100 µS/cm, 5x60 mL bottles
011007	Orion conductivity standard, 1413 µS/cm, 5x60 mL bottles
011006	Orion conductivity standard, 12.9 mS/cm, 5x60 mL bottles
1010001	Orion conductivity calibration resistor kit

For more information, contact your local Thermo Scientific products dealer or call our customer and technical service experts at 1-800-225-1480 (for the US and Canada) or visit www.thermoscientific.com/water.

©2011 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. ROSS and the COIL traddress are trademarks of Thermo Fisher Scientific Inc. US Patent 6,793,787.



S-STARA215-E-1011-RevB

Water Analysis Instruments
www.thermoscientific.com/water

North America
166 Cummings Center
Beverly, MA 01915 USA
Toll Free: 1-800-225-1480
Tel: 1-978-232-6000
info.water@thermo.com

Netherlands
Tel: (31) 033-2463887
info.water.uk@thermo.com

China
Tel: (86) 21-68654588
wai.asia@thermofisher.com

India
Tel: (91) 22-4175-8800
wai.asia@thermofisher.com

Singapore
Tel: (65) 6778-6876
wai.asia@thermofisher.com

Japan
Tel: (81) 045-453-9175
wai.asia@thermofisher.com

Australia
Tel: (613) 9757-4300
infowaterau@thermofisher.com

Thermo
SCIENTIFIC